

Wall Mounted Workstation

User Manual

For Non-Locking and PIN Code Locking Cabinet Models



Warnings



IMPORTANT – Indicates a situation that does not present any hazard but is very important in maintaining a well functioning workstation.



ATTENTION – Consult manual to avoid a potentially hazardous situation which may result in minor or moderate injury.



ELECTRICAL – Indicates an impending electrical hazard which, if not avoided, may result in personal injury, fire and/or death.




- Contact the Facility Engineer for direction on mounting locations and methods prior to installing any wall tracks or equipment
- The shipping weight of this unit is 86 lbs (39 kg). Use proper lifting techniques to prevent injury
- The supplied power cord is rated for medical use. Connecting the cord to an outlet that is not medical grade (indicated with a green dot) will not ensure grounding protection (locking cabinet only).
- Power cord, USB extension, and workstation are for INDOOR use only. **DO NOT OPERATE OUTDOORS.**
- Keep power cord away from water. **DO NOT PLUG CORD INTO OUTLET IF WET.**
- **DO NOT OPERATE PRODUCT IF WET.** If the WORKSTATION becomes wet, unplug it immediately, wipe off any excess liquid, and allow it to dry before using again.
- Inspect power cord before integration. **DO NOT USE POWER CORD IF DAMAGED.**
- Fully insert power cord plug into outlet. **DO NOT** unplug by pulling on cord. **DO NOT** remove, bend or modify any metal prongs or pins of power cord.
- **DO NOT** use excessive force to make mechanical or electrical connections
- **DO NOT** use an electrical extension cord with your workstation
- **DO NOT** use a flammable cleaner on the station as it can result in fire or explosion

Table of Contents

Warnings	2
Installation	4
Box Contents	4
Parts Breakdown.....	5
Wall Cabinet Disassembly.....	6
Mounting Wall Cabinet.....	7
Wall Tracking Installation.....	8
Wall Cabinet Reassembly	8
Integration	9
Technology Schematic.....	9
Specifications	10
Technology Integration.....	11
USB Cable Locations	13
Software Installation	14
USB Task Light Software	14
Programmable Controller/Lock	16
Locking Cabinet PIN Code	16
Maintenance	18
Cleaning	18
Wood Panel Care	18
Troubleshooting	19
Warranty	20
Limited Warranty for Wall Mounted Workstations	20
Service Details.....	20
Transport/Storage	21
Statement of Use	21
Service	21
Service Request	21
Service Level Commitment	21
Standards Compliance	22


Installation

BOX CONTENTS

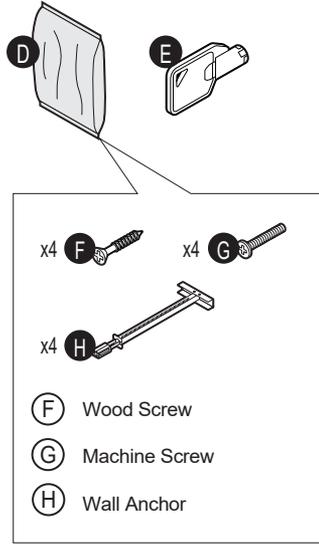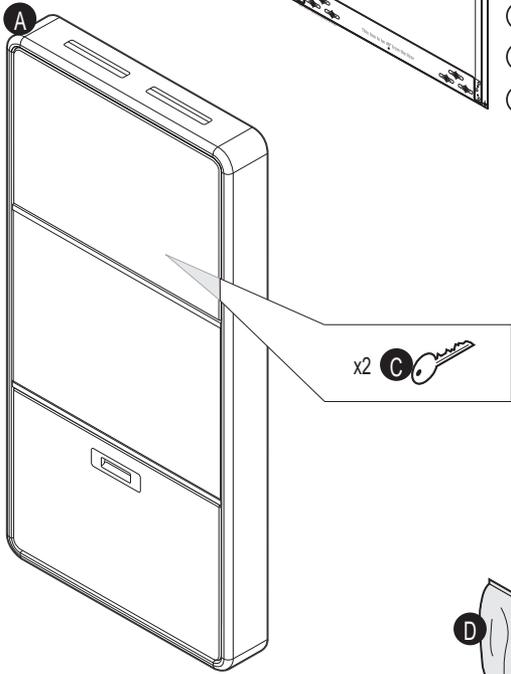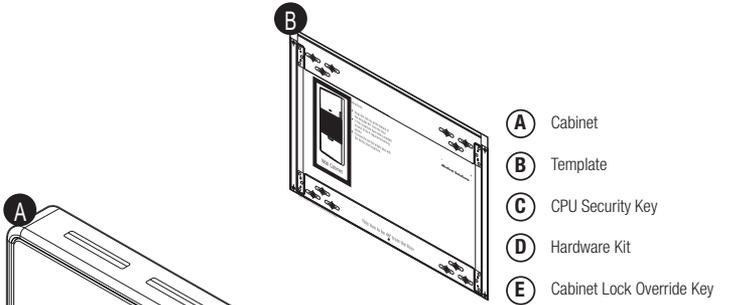





Installation

PARTS BREAKDOWN

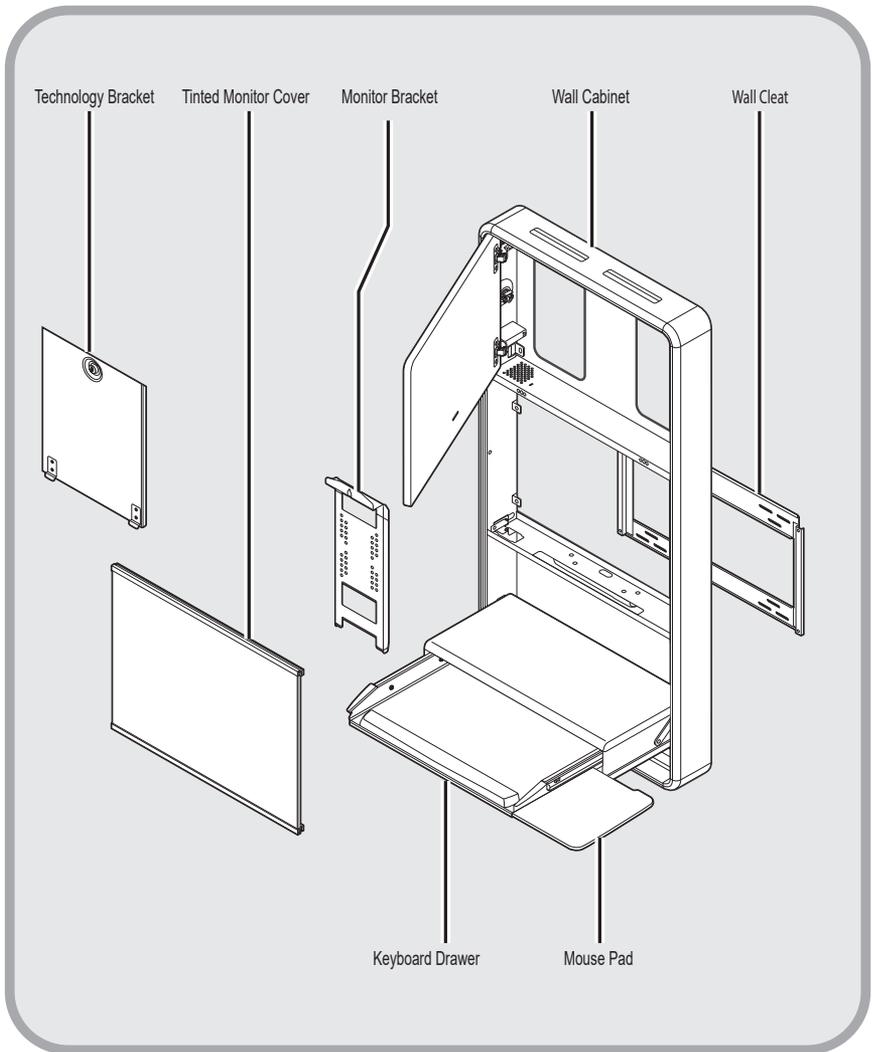


Installation

WALL CABINET DISASSEMBLY

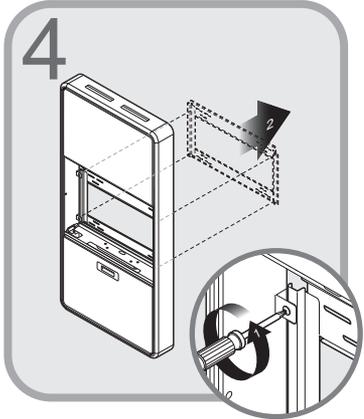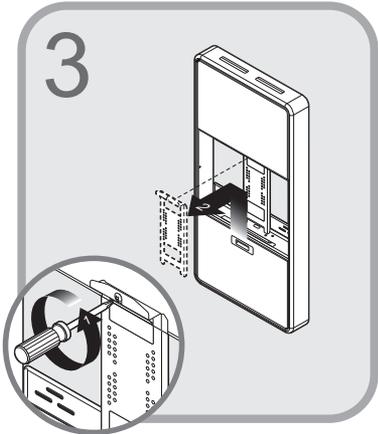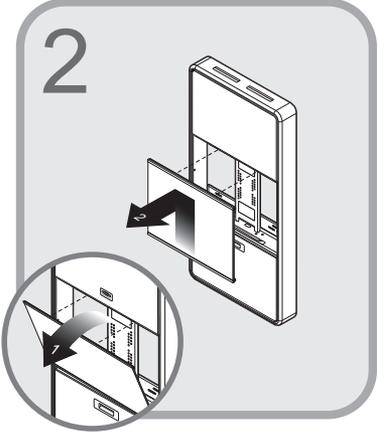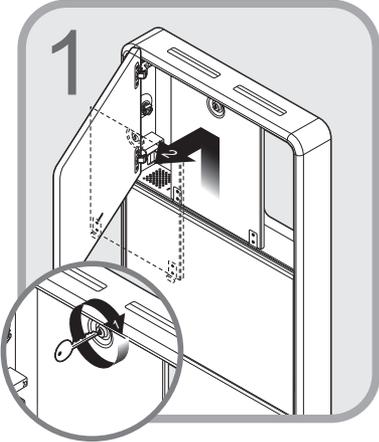


Installation

MOUNTING WALL CABINET

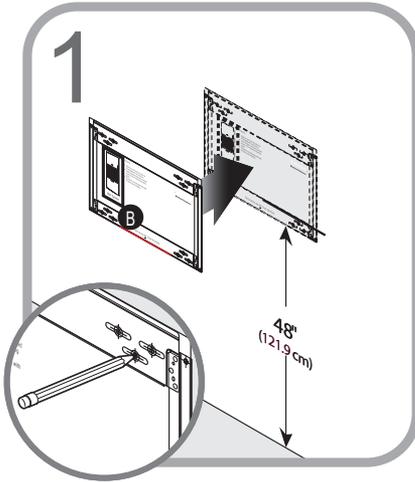

Note: Consult a Facilities Engineer regarding structural codes and utilities.

Use the mounting template (B) to determine drill hole locations. Set line should be 48" (121.9 cm) from floor. Mark and drill holes.

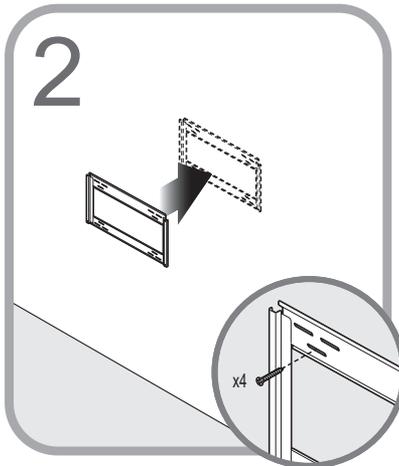


Mount wall cleat using hardware provided.

Note: Use part wall anchor (H) and machine screw (G) for steel stud installation. Use wood screws (F) for wood studs.



Installation

WALL TRACK INSTALLATION

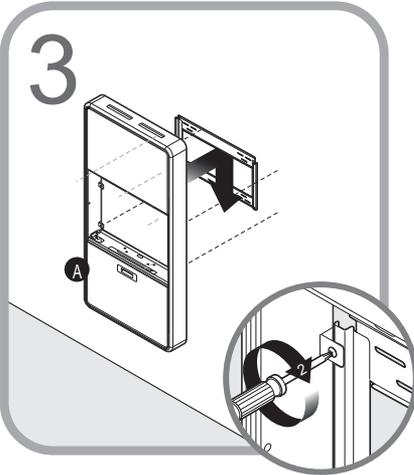


Hang wall cabinet (A) on wall cleat and secure with screws.

WALL CABINET REASSEMBLY

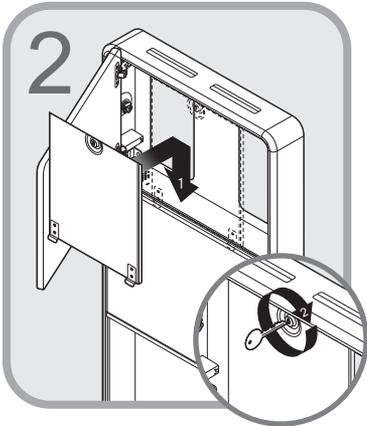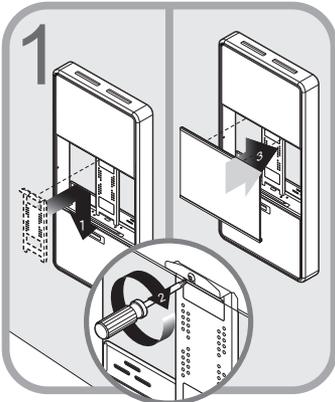

Integration

TECHNOLOGY SCHEMATIC

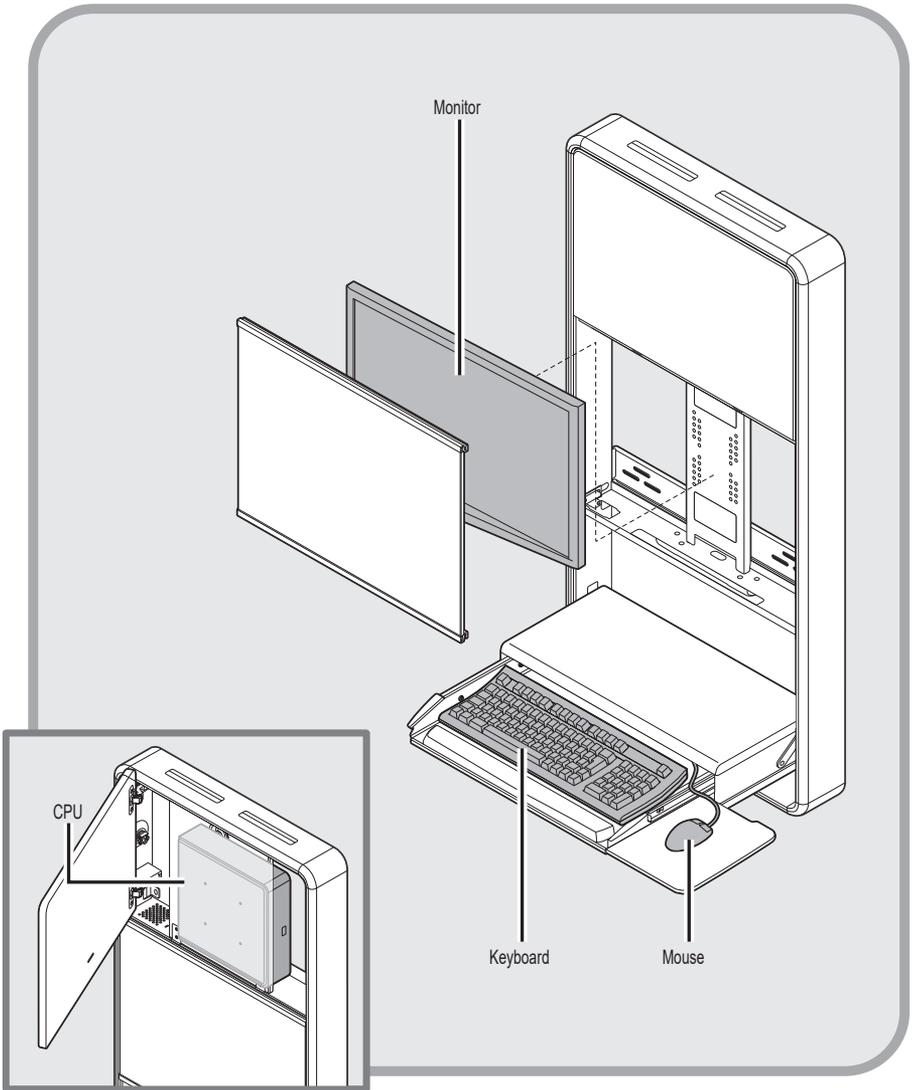





Integration

SPECIFICATIONS

Power Cord (locking cabinet only)	1.5' (45.0 cm) length, medical grade NEMA right angle plug, 3x NEMA receptacles, 125 V 50/60 Hz, 250 W max
USB Cables	2 – USB 2.0 type A male to type A female extension cables
Monitor Bracket	25 lbs (11.3 kg) max; VESA 75 mm & 100 mm
Keyboard Drawer	Accommodates 1.75" H x 18" W x 8" D (4.5 cm x 45.7 cm x 20.3 cm) USB keyboard
Work Surface	20.38" W x 10" D (51.8 cm x 25.4 cm)
Mouse Pad	10" W x 10" D (25.4 cm x 25.4 cm)
Technology Cabinet	Accommodates CPUs up to 11" W x 11" H x 2.9" D (27.9 cm x 27.9 cm x 7.4 cm)





Integration

TECHNOLOGY INTEGRATION

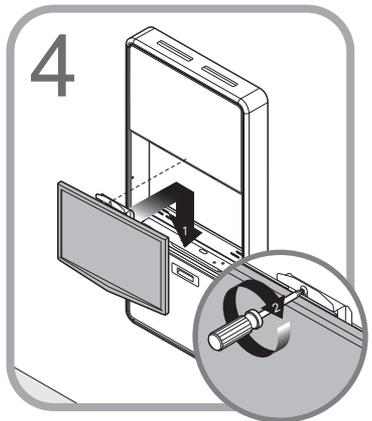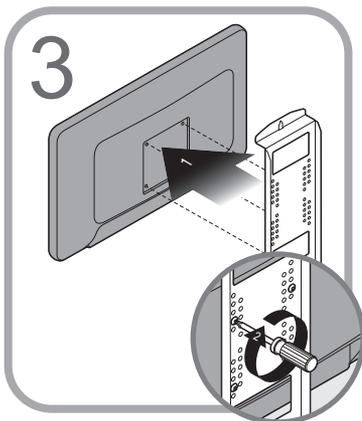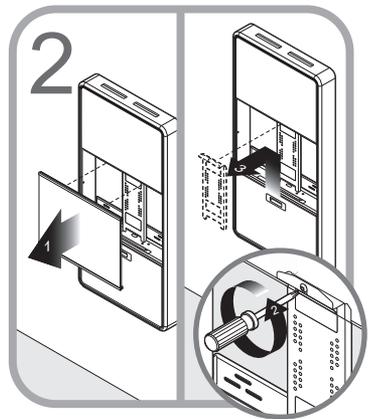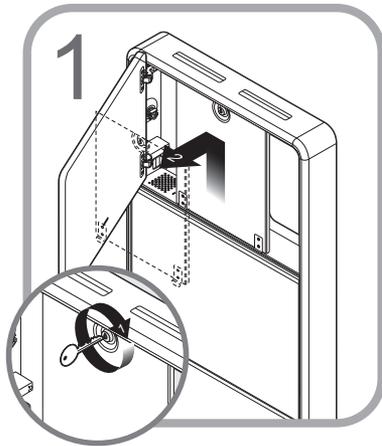


Integration

TECHNOLOGY INTEGRATION

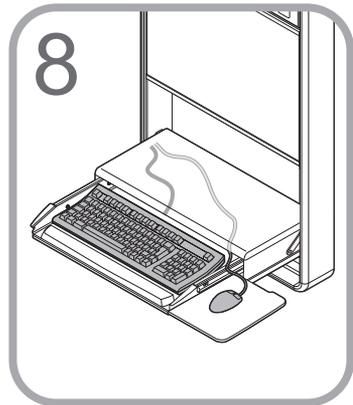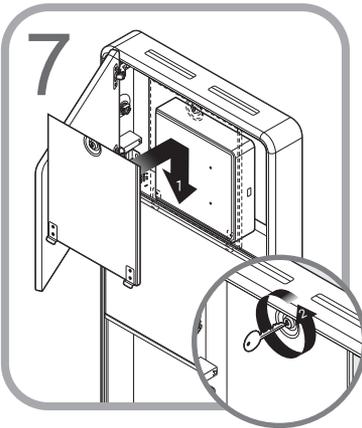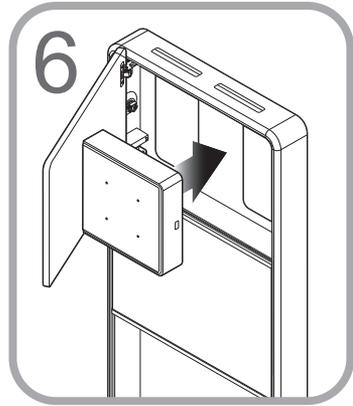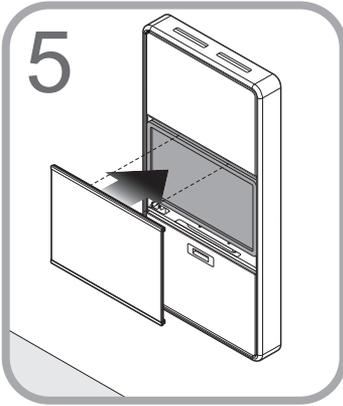

Integration

USB CABLE LOCATIONS

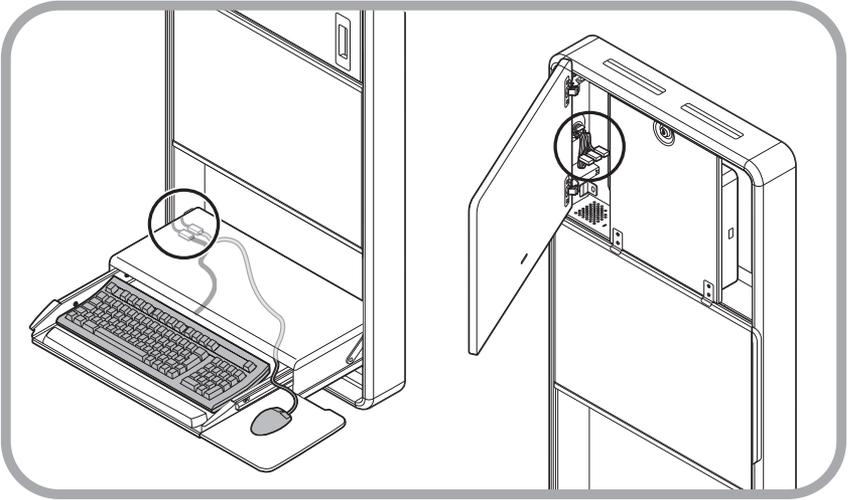

Two USB cables are located in the rear corner of the keyboard tray. These are routed to the technology box area and should be plugged into the CPU.

Three USB cables are located in the technology box area. These should be connected to CPU to provide power for the Task Light, Keyboard, and Monitor.

Software Installation


USB TASK LIGHT SOFTWARE

Preparation


Uninstall any previous drivers that may have been installed for the USB task light.

Note: If the USB device has not been installed previously, skip to *Installation* section on page 15.


1. Start > Control Panel > System.
2. Click on the *Hardware* tab.
3. Select the **Device Manager** button.

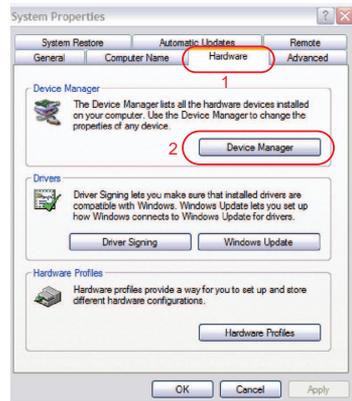

4. Scroll down to *Ports (Com & LTP)*.
5. Expand the dropdown and select *AT90USBxxx CDC USB to UART MGM (Com XX)*.

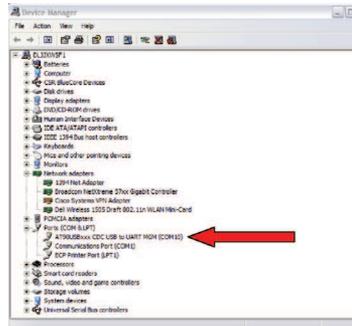

5. Click the *Driver* tab at the top (1).
6. Select *Uninstall* (2).
7. Click on **OK** (3).
8. Restart the computer.

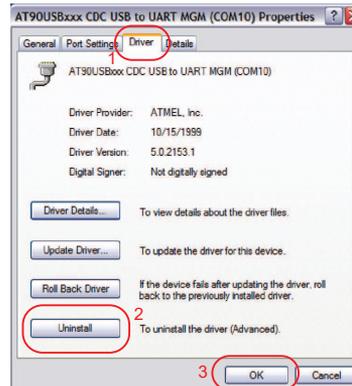

Note: If the drivers did not install correctly the first time it may be listed in the *Other* device list usually indicated with a yellow Asterisk * symbol.

Software Installation

USB TASK LIGHT SOFTWARE

Installation

Note: DO NOT plug the USB cable in at this point.

1. Run the *RubbermaidLCDSecureIT.exe* setup file.
2. Once installation is complete, plug in the USB to an open port and wait for the driver request.

Note: The USB driver is located within the files provided.


3. Select *at90usbxxx_cdc.inf*.
4. After the USB drivers have been installed, you can setup the software preferences (see below).

Software Setup

There are 4 options that are available when the keyboard is closed:

1. **None** (No action take)
2. **Immediate Blank Screen**
3. **Immediate Screen Saver** (select *Resume*, *Password Protect* under *Screen Saver Options* if user ID and password is required)
4. **Immediate Lock Computer**

Note: You can adjust the duration of time the light (5) will stay on at any given instance. The default is set at 5 minutes. The password section is not applicable in this situation (non-locking cabinets).

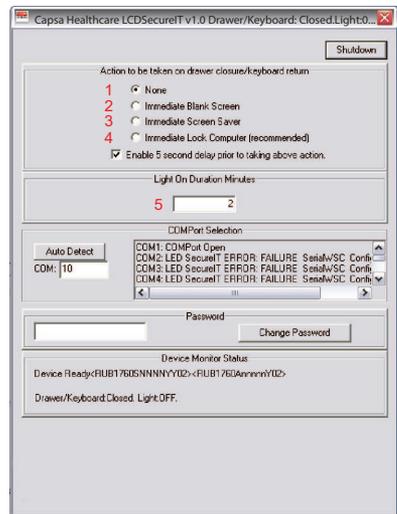


Shutdown the software

In the upper right-hand corner there is a shutdown  button to completely close the software.

Note: If you close the window, it will be minimized to the system tray.

Programmable Controller/Lock

LOCKING CABINET PIN CODE

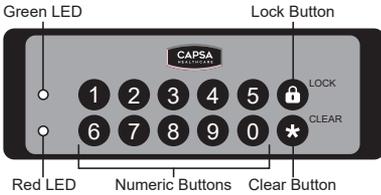

In normal operation the user will enter a 4-digit access code. Once a numeric key is pressed, the green LED will light. The user will then enter the remaining 3-digits of their code. If the access code is accepted, the green LED will blink and the door actuator will be released for the time period specified. If the password is not accepted, the red LED will blink.

If a mistake is made during code entry, the '*' key can be used to clear previous key strokes and start over. This will also turn off the green LED. Once a key is pressed, the keypad will wait 30 seconds for the next input. If 30 seconds elapse with no input, the keypad will clear previous key strokes. The green LED will turn off and the red LED will flash.

Note: Once a function code or supervisor code has been accepted, the keypad will wait 30 seconds for the next key input. If 30 seconds elapse with no input, the red LED will flash and the keypad will return to normal mode. The **Lock** key will be inactive during keypad programming mode. Serial communication will be inactive during keypad programming mode. During serial communication, all keypad buttons will remain inactive.

Programming Functions

All programming functions are initiated by entering the * key, then the supervisor code, followed by the * key. Next, the programming code is entered, followed by the * key. See the table below for programming codes.

Programming Codes

0*	Restore default settings
1*	Erase all user codes
2*	Remove specified user code
3*	Add user code
4*	Program re-lock time
5*	Program new supervisor code
6*	Program keypad serial number

LED Light Status

Solid Green	Waiting for user input
Flashing Green	Released door, command accepted
Solid Red	Serial communication, keypad inactive
Flashing Red	Lock command, input error

Unlocking the Control Pad



The green LED will flash and the door actuator will be released.

Restore Default Settings

The controller can be set to its factory default settings. To restore the default settings:



The green LED will flash indicating that programming is complete. The table below lists the default parameters.

Default Parameters

Default Settings	
User code	8294
Supervisor code	7183
Door re-lock time	10 seconds
Keypad serial number	00-00000-WC-01

Note: Restoring the default settings will also erase all user codes except the default code.

Programmable Controller/Lock

LOCKING CABINET PIN CODE

Erase All User Access Codes



The green LED will flash indicating that erasing is complete.

Remove User Access Code



The green LED will flash indicating that erasing is complete. If the code is not stored in memory, or an invalid key is pressed (*, 'Lock'), the red LED will flash.

Add User Access Code

128 user access codes can be programmed into the keypad. The user access codes can be any combination of digits. To program an access code:



The green LED will flash indicating that programming is complete. If the user memory is full, or an invalid key is pressed (*, 'Lock'), the red LED will flash. If the code entered is already stored in memory, the green LED will flash as normal.

Program Keypad and Serial Number



The green LED will flash indicating that programming is complete. If an invalid key is pressed (*, 'Lock'), the red LED will flash.

Program Door Release Time

The length of time the door is released can be programmed from one second to 90 minutes. To change from the default of ten seconds:



The green LED will flash indicating programming is complete. If an invalid key is pressed (*, *, 'Lock'), the red LED will flash.

Note: Enter the desired amount of time in the format 'MMSS' where :

MM is the number of minutes and **SS** is the number of seconds.

Maximum time allowed is 90 minutes. Any time above the maximum will default to 90 minutes.

Program Supervisor Code



The green LED will flash indicating that programming is complete. If an invalid key is pressed (*, 'Lock'), the red LED will flash. If the code entered is already stored in memory, the green LED will flash as normal.

Maintenance



DO NOT use the workstation if pieces are missing or the unit is damaged. In these cases, immediately contact Capsa Healthcare Customer Service for more information: 800-437-6633.

Cables - Always keep the cables neatly organized and be sure to route cables away from moving components with wire ties or cable clips.

Power cord (locking cabinet only) - Periodically inspect power cord and plug to ensure plug is not bent and cable is not frayed.

CLEANING



CAUTION: Because of the close proximity of electrical power and equipment, flammable cleaners should never be used on the workstation.

- Verify that your workstation is unplugged from the wall outlet before cleaning
- Allow your workstation to dry completely before plugging the power cord into a wall outlet
- When cleaning the workstation, wipe surface with a damp cloth and thoroughly dry.
- Never cover the workstation or its components with liquid or allow liquids to flow into the workstation
- Never use steel wool or other abrasive material as these could damage the surface finish
- Before using any cleaner on the workstation, first test on a small area to ensure that the surface is not harmed
- These guidelines cannot guarantee infection control. The hospital's Infection Control Administrator should be consulted regarding cleaning procedures and schedules.
- Clean plastic components with diluted, non-abrasive solutions. Suggested cleaners are water, soap, diluted bleach and alcohol solutions.
- Remove pen and dry erase marker stains with a soft cloth and 91% isopropyl alcohol
- Remove iodine stains with a soft cloth and any cleaners suggested above



DO NOT use the following chemicals to clean your workstation: acetone, mineral spirits, abrasive cleansers, paint thinner or any other harsh or toxic chemicals.

WOOD PANEL CARE

Harsh, abrasive and undiluted cleaning products may cause damage to Deco Lam® and the contact adhesive.

Pine Sol® and Simple Green® are cleaners that have been approved for use on Deco Lam®. A 30-1 ratio of water to cleaner is highly advised. Recommend water and a clean towel as a cleaning alternative to homeowners.

Troubleshooting

Problem	Solution
Keyboard mechanism folds down too fast or too slow.	Remove the right side extrusion and locate the adjustment screw on the damper. To slow it down, turn screw to right. To make it faster, turn screw to the left.
Top door needs adjustment for alignment.	Open the top door and locate the 2 hinges on the left side. To adjust the door up or down use the inside screw of the hinges. To adjust the door right or left use the outside screw on the hinges.
Task Light is not turning on.	Remove the right side extrusion and make sure the light assembly cables with orange ends are secure on the switch.

Warranty

LIMITED WARRANTY FOR WALL MOUNTED WORKSTATIONS

Capsa Healthcare is pleased to offer a five-year warranty on durable components and a two-year warranty on electronic components.

If during the warranty period this Capsa Healthcare product proves defective in materials or workmanship under normal use by the original purchaser, please contact Capsa Healthcare technical support at www.capsahealthcare.com/service (please be sure to complete all information, including product serial number, description of the issue, and full contact information). Capsa Healthcare will determine, at its sole discretion, how to best address your warranty issue, which may include sending you a replacement part covered under warranty or for sale. Capsa Healthcare reserves the right to require proof-of-purchase prior to honoring any warranty request. This warranty does not cover product abuse, modification, failure to adhere to product instructions, or improper operation/misuse. Capsa Healthcare SHALL NOT BE LIABLE FOR ANY CONSEQUENTIAL OR INCIDENTAL DAMAGES WHATSOEVER. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you. This warranty gives you specific legal rights and you may also have other rights which vary from state to state or country to country.

SERVICE DETAILS

Consumable components are not covered under warranty and include:

- Side Extrusions
- Cable anchors
- Backboard
- Mouse Tray

All other standard components will be replaced under the applicable warranty following a filed service request.

Transport/Storage

- The shipping weight of this unit is approximately 86 lbs (39 kg). Use proper lifting techniques to prevent injury.
- Care should be taken to transport and store this system within a temperature range of 32° F to 90° F (0° C to 32° C); humidity 20% RH to 95% RH non-condensing.

STATEMENT OF USE

The Wall Mounted Work Station is designed and manufactured by Capsa Healthcare. Capsa Healthcare's goal is to elevate your performance by offering innovative solutions that positively impact clinical effectiveness and enhance medication management processes throughout all healthcare environments.

- The Wall Mounted Work Station is designed for safe use in general patient areas for the purpose of clinical data entry and retrieval
- The Wall Mounted Work Station is not intended for home use
- The Wall Mounted Work Station has no potential electromagnetic or other interference risks when operated according to guidelines covered in this instruction manual

Service

SERVICE REQUEST

Contact your IT department, or file a service request at our website:
www.capsahealthcare.com/service.

SERVICE LEVEL COMMITMENT

Capsa Healthcare is committed to providing best-in-class service. This document details our standard warranty and instructions on how to request service using our customer support system.

Standards Compliance

TEST COMPLIANCE

NRTL

- o UL 60601-1 - Issued:2003/04/25 Ed:1 Rev: 2006/04/26 UL Standard for Safety Medical Electrical Equipment, Part 1: General Requirements for Safety
- o CSA C22.2 No. 60601-1 - Issue: 2008/02/01 Ed:2 Medical electrical equipment - Part 1: General requirements for basic safety and essential performance; COR 2: 2011/06/01
- o IEC 60601-1-1 - Issued:2000/12/01 Ed:2 Medical Electrical Equipment - Part 1-1: General Requirements for Safety - Collateral Standard: Safety Requirements for Medical Electrical Systems

IEC/EN 60601-1-2 (Ed. 2): 2001 +A1: 2004 - Medical electrical equipment Part 1-2: General requirements for safety - Collateral standard: Electromagnetic compatibility

- o Conducted Emissions - Conducted Voltages (CISPR 11/ EN 55011 (CV))
- o Radiated Emissions - Electric Fields (CISPR 11/ EN 55011 (RE-E))
- o Electrostatic Discharge Immunity Test - (IEC 61000-4-2)
- o Radiated, radio-frequency, electromagnetic field immunity test (IEC 61000-4-3)
- o Electrical Fast Transient/Burst Immunity Test (IEC 61000-4-4)
- o Surge Immunity Test (IEC 61000-4-5)
- o Immunity to conducted disturbances, induced by radio-frequency fields (IEC 61000-4-6)
- o Power Frequency Magnetic Field Immunity Test (IEC 61000-4-8)
- o Voltage Dips, Short Interruptions and Voltage Variations Immunity Tests (IEC 61000-4-11)

FCC

PART 15, Subpart B, Class A – Unintentional Radiators

PRODUCT CLASSIFICATION

- o Class 1/ Internally powered device with no applied parts.
- o This equipment is designed for continuous operation.
- o Class A, Group 1 ISM Equipment
- o This device is classified IPX0 for water ingress
- o Input 500 Watts North America
- o Input Voltage 120 Vac 60 Hz

Notes



800.437.6633
www.CapsaHealthcare.com