


Automated Narcotics Storage, Access, & Tracking



Ask Administration or Nursing About NexsysADC™

- Q1:** How do you store controlled substances? In a safe or closet?
- Q2:** Is drug diversion a concern?
- Q3:** If you considered Automated Dispensing Cabinets (ADCs), were they too expensive or complex?
- Q4:** Are you prepared for a DEA audit? Do you have electronic records?

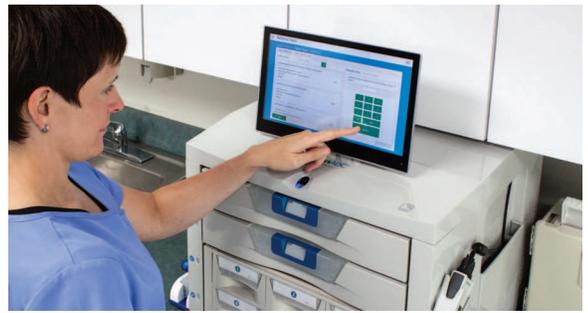

Busy Ambulatory or Oral Surgery Centers, Acute Care and Other Care Sites

When medications are stored in NexsysADC, you always know the inventory, who is accessing the cabinet, and what is dispensed for a specific patient or procedure. NexsysADC is the practical solution to:

- Increase visibility to inventory onhand
- Deter internal diversion, especially controls
- Ensure medications are onsite when needed
- Promote efficient patient care

The Practical ADC For Any Care Facility Storing Controlled Meds

NexsysADC, from Capsa Healthcare, is the simplest system to closely manage the onsite storage of controlled or high-value medications and supplies.

- Costs 50% to 75% less than other automated dispensing cabinets
- Limits & tracks access to narcotics and high-value meds and supplies
- Main cabinet and tabletop sizes
- Select from 4 popular configurations -- See p. 3



4
New
Packages
(Page 3)



4T Countertop Cabinet
For smaller spaces,
budgets, formularies



Main Cabinet
For busy facilities with a
higher patient population

Contact Capsa Healthcare for more information
lfc.support@capsahealthcare.com
800.437.6633

New!

4 NexsysADC Popular Configurations

(see page 3)

1



2



3



4



Automated Narcotics Storage, Access, & Tracking



How NexsysADC™ Works

NexsysADC™ is available as a free-standing 10-tier main cabinet or a countertop size 4-tier cabinet. In either size, the ADC will store, track, and dispense controlled medications, and high-value meds and supplies. Staff member access is controlled with biometric, or proximity card, or typed ID. With intuitive software and pick-to-light technology, the clinician finds the exact medication or supply needed, and NexsysADC automatically logs the transaction and updates the onboard inventory.

Each NexsysADC has 2 or 3 Storage Zones for controlled and high-value medication and supplies



1 Controlled Medications

Controlled and high-value medications are safely guarded in CAMs (controlled access modules) with tamper-evident locking lids. Choice of 3 CAM layouts.

2 Open Cell Medication Bins

Routine meds and high-value meds are efficiently managed in dividable, locked open-matrix bins, with pick-to-light guidance. Each bin configurable for 1-14 items.

3 Bulk Supplies & Kits

Bulk supplies and larger items are stowed neatly in configurable, locked supply drawers, with pick-to-light guidance. Each drawer configurable for 1-16 items.

NexsysADC™ CAMs & Features


Unique CAMs are at the heart of NexsysADC to protect controlled medications. Illuminated prompts direct you to the correct dose. Includes clear locking lids and tamper evidence.



Depending on your configuration (p. 3), CAMs hold 11, 13, or 19 SKUs of any dosage form (oral dose, syringe, ampule, vial) in individual cells.



NexsysADC is simple to navigate, with clean, bright software design and intuitive user interface.

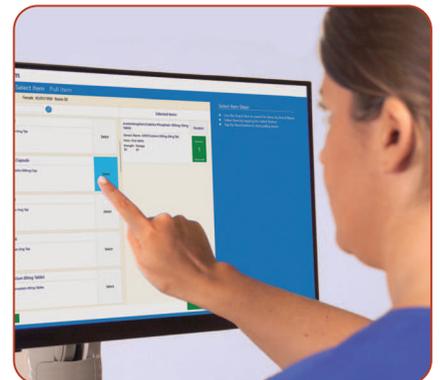

Contact Capsa Healthcare for more information
lrc.support@capsahealthcare.com
800.437.6633

New! 4 NexsysADC Configurations

If these configurations don't suit your needs, contact Capsa for a custom layout



1

Configuration 1, Main Cabinet

For Controlled Meds: CAM-11 holds 11 SKUs

For Non-Control Meds & Supplies: Two 3-tier dividable open-matrix inserts (18 bins); one 3" and one 6" dividable, open-matrix supply drawers (self-configurable, total storage up to 284 SKUs)

Part # NX-MAINBDL-1



2

Configuration 2, Main Cabinet

For Controlled Meds: CAM-11 holds 11 SKUs; CAM-13 holds 13 SKUs

For Non-Control Meds & Supplies: Two 3-tier dividable open-matrix inserts (18 bins); one 6" dividable, open-matrix supply drawer (self-configurable, total storage up to 268 SKUs)

Part # NX-MAINBDL-2



3

Configuration 3, 4T Countertop Cabinet

For Controlled Meds: CAM-11 holds 11 SKUs

For Non-Control Meds & Supplies: One 3-tier dividable open-matrix insert (9 bins, self-configurable, total storage up to 126 SKUs)

Part # NX-4TBDL-1




4

Configuration 4, 4T Countertop Cabinet

For Controlled Meds: CAM-19 holds 19 SKUs; CAM-11 holds 11 SKUs

For Non-Control Meds & Supplies: One 2-tier dividable open-matrix insert (6 bins, self-configurable, total storage up to 84 SKUs)

Part # NX-4TBDL-2

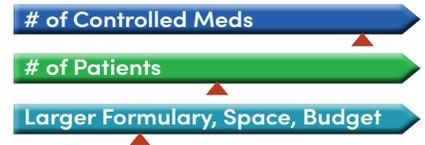

CAM-11
11 cells
Included with configuration 1, 2, 3, 4



CAM-13
13 cells
Included with configuration 2



CAM-19
19 cells
Included with configuration 4

Specifications & Components

4T Cabinet Size

23.5"(w) 20.4"(d) 17.7"(h) with 14" touchscreen monitor. Approximately 100 lbs. empty. 4 Tiers

Main Cabinet Size

25"(w) 20.4"(d) 43"(h) with 22" touchscreen monitor. 145 lbs empty. 10 Tiers

Includes:

Biometric scanner, keyboard, barcode scanner for medication verification, interfacing capability, implementation & training

Contact Capsa Healthcare for more information

ltc.support@capsahealthcare.com
800.437.6633