

Avalo® Series Medical Cart Line

Emergency Preparedness Cart

Today, emergency preparedness is more important than ever. Infectious disease control requires its own, unique approach to workflow processes, supplies organization, and medication control. The Avalo Emergency Prep Cart design includes features that promote simple cleaning and post-use infection control:

- Smooth panels and work surfaces
- Durable hi-impact panels & integrated bumpers
- Seamless, molded drawers
- Simple to remove / replace drawer dividers


Flexible drawer configurations and a wide selection of accessories enhance user efficiency and support emergency preparedness protocols.




Accessories & Options



I.V. Pole

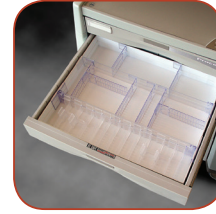

Removable Tray

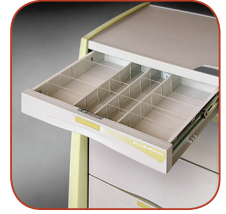

Drawer Dividers



Storage Module



Waste Container



Keyless Locking

[Contact Capsa to inquire about additional Accessories & Options](#)

Cart Specifications

	Surface Height	Depth	Width
Standard	43.5" / 111 cm	24" / 61 cm	31" / 79 cm
Intermediate	39.5" / 100 cm	24" / 61 cm	31" / 79 cm
Compact	36" / 91 cm	24" / 61 cm	31" / 79 cm

Drawer Specifications

	Height	Depth	Width
Three Inch	3" / 7.6 cm	17.5" / 45 cm	21.5" / 55 cm
Six Inch	6.25" / 16 cm	17.5" / 45 cm	21.5" / 55 cm
Ten Inch	10" / 25 cm	17.5" / 45 cm	21.5" / 55 cm

Model Options (custom configurations available)



Standard



Intermediate



Compact

Specifications, availability & components are subject to change without notice.

Emergency Preparedness Cart Features

- Precision-welded steel frame with durable, high-impact panels
- Lightweight design and profile simplify movement and steering
- Seamless drawer trays are easy-to-clean and feature premium, self-closing drawer slides
- Wide selection of accessories and organization systems
- Expansive work area with optional slide-out work surface
- Delivered fully assembled, no tools required
- 10-year limited warranty